

PRISM

Prism Hammock Sling Range



User Manual

PRISM MEDICAL UK

Prism Hammock Sling

Application

The Hammock sling provides maximum support to the sacral area and is used in cases where the Universal Sling is unsuitable because the individual has poor lower body tone. It can be used to transfer from the seated to recumbent position or from recumbent to seated. The person must be assessed by a competent person who must be fully trained in the suitability, application and fitting of the sling. Carers should always be trained in the use of hoists and slings.

This sling is usually used in a divided leg position however can be used in a cross over position for amputees or other people following an accurate risk assessment.

Sling Inspection

All Prism slings are manufactured to meet and exceed the requirements of BS EN 10535:2006. Slings should be inspected prior to use and after washing. Damaged or badly worn slings should be discarded.. The label on the sling contains vital information to identify the sling - if any part of the label becomes illegible, then it must be removed from service and replaced.

The sling should be visually inspected with regard to the following items:-

- Any signs of damage (i.e. Cuts, frays, tears, burns) to the sling straps.
- Any signs of stitching becoming undone on any part of the sling.
- Any signs of the fabric being stretched, tearing or undue wear.
- Any signs of discolouring or bleaching on the material.
- Is the sling complete? i.e. All straps present, stiffeners in place - where applicable, etc.
- Label clearly legible.

Compatibility

The issue of sling compatibility with various manufacturers' hoists has been a concern for specifiers for a number of years, an issue that is compounded by the increased number of companies supplying hoists and slings.

Some hoist/sling manufacturers will state that only slings manufactured by them can be used on their hoists and vice versa. This is their policy only and is **NOT** a legal requirement. A risk assessment must still be completed even if the same manufacturer of hoist and sling are to be used, and the specifier simply has to state that the carry bar is of a similar design and that the sling manufacturer has authorised its use on it.


A comprehensive list is available separately and your local agent will be able to provide assistance with compatibility issues.

Maximum Load (Safe Working Load - S.W.L)

The maximum load for the sling is clearly marked on the label of the sling. The maximum load will also be displayed on the hoist and any detachable elements of the hoisting system.

DO NOT exceed the maximum load stated on any of the equipment.

Sling Dimensions

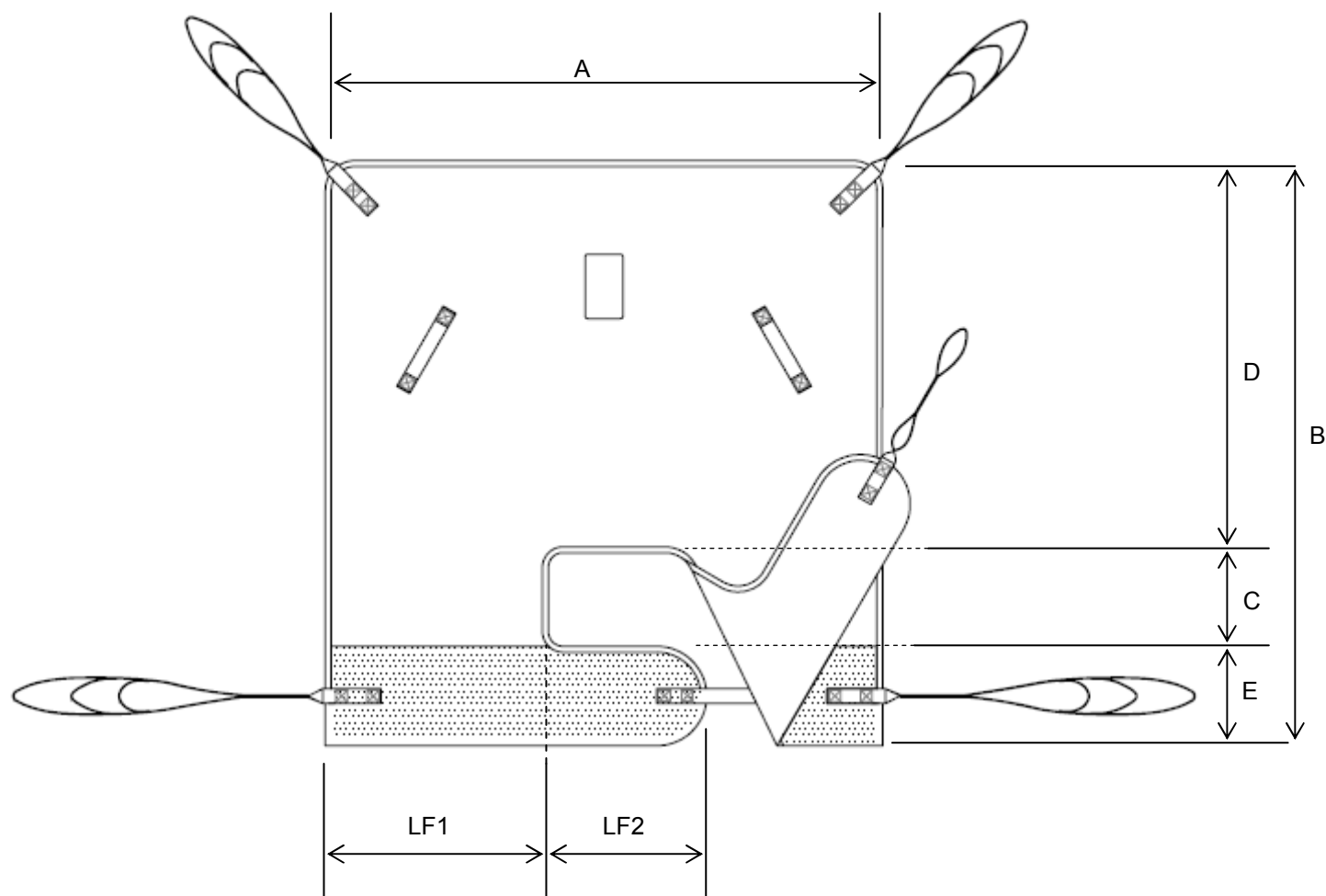


	A	B	C	D	E	LF1	LF2
Child	640	670	115	450	105	245	200
Junior	750	850	150	580	120	310	240
Small	850	950	160	630	160	310	350
Medium	960	1060	185	700	175	370	360
Large	1050	110	200	700	200	390	430
X Large	1180	1180	250	730	200	400	495


All dimensions have a tolerance of +/- 20mm

Hoisting from a chair


Note: The person to be transferred shall be referred to as “the person”.

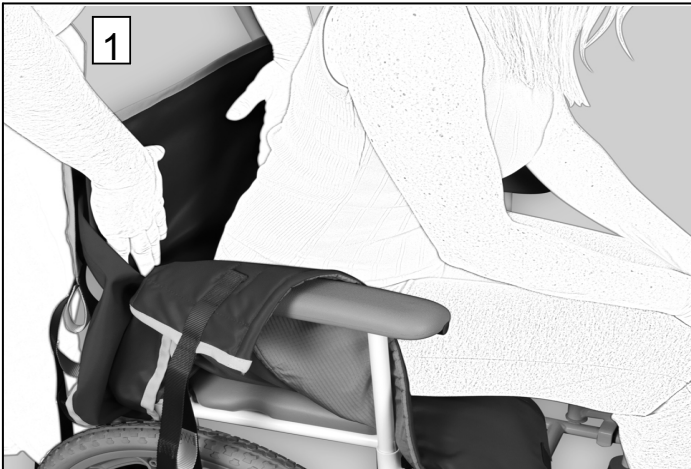


Introduce the sling by sliding it down between the chair backrest and the person's back until the top of the commode aperture reaches the base of the spine. If the sling can not be introduced from behind, apply the sling in a position level with the person. Leave the sling half way down the person's back (at this stage don't bring the sling to shoulder level).

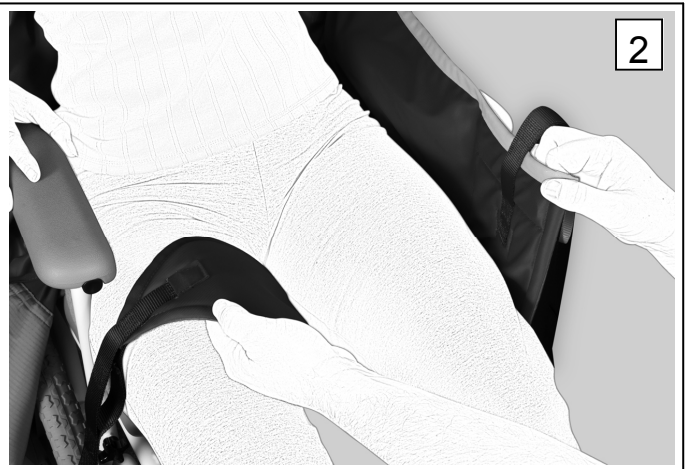

Bring the leg sections along-side the person's thighs and pull lightly to ensure that they are of equal length. Using good manual handling techniques, raise the person's leg and feed the leg sections under the person's thighs.

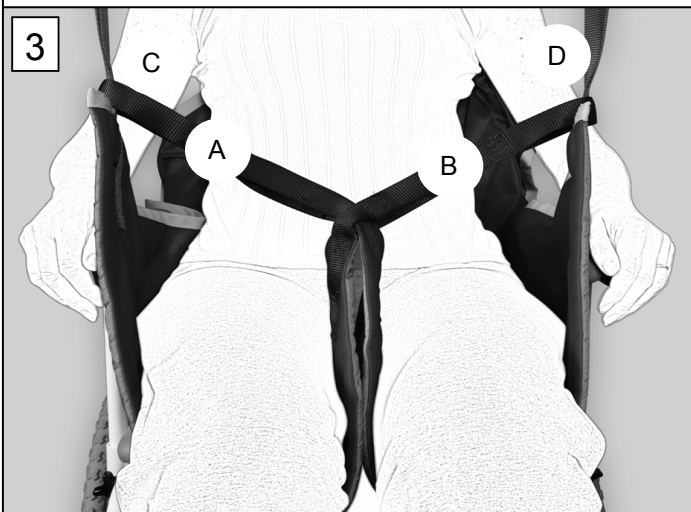

Cross strap A through the loop at the base of strap B. Pass strap C through loop A and pass strap D through loop B. Now raise the top of the sling level with the person's shoulders. The fitting of the sling is now complete.

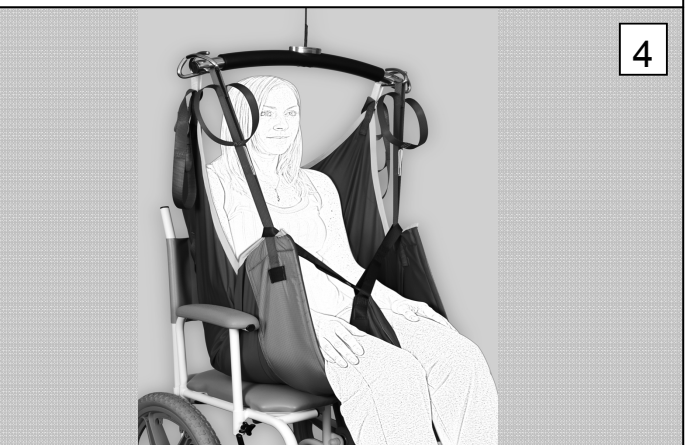


The four straps should now be attached to the carry bar hooks. Raise the carry bar just enough to tension the straps and **ensure that they are still securely attached**. The person can now be lifted as required. To remove the sling, carry out the fitting procedure in reverse. Important - always peel the sling away from the person - this avoids skin friction and movement.

To remove the sling, carry out the fitting procedure in reverse.
Important - always peel the sling away from the person to avoid skin friction and movement.

Hoisting from a lying position




Using good manual handling techniques roll the person towards you so that they are positioned on their side. To introduce the sling, roll half of the sling up underneath itself so that when the person is rolled onto their opposite side, the sling can be unrolled until it is flat. Use the centre markers to ensure the base of the back section is in line with the base of the person's spine. Ensure that handle loops are to the outside.

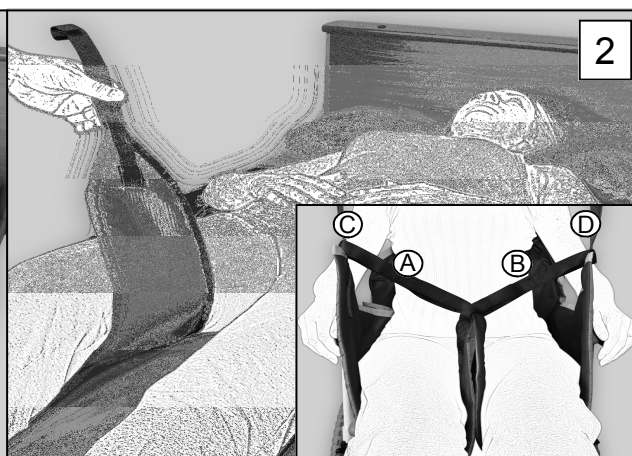

Pass the leg section under and up between the thighs. Ensure that as much section as possible is pulled under the leg and that all creases are eliminated to achieve maximum comfort.

Cross strap A through the loop at the base of strap B. Pass strap C up through strap A and pass strap D through strap B. The fitting of the sling is now complete.

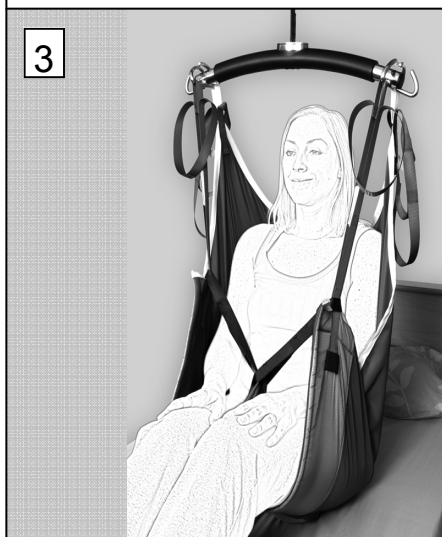

The four straps should now be attached to the carry bar hooks. Raise the carry bar just enough to tension the straps and **ensure that they are still securely attached**. The person can now be lifted as required.



Hoisting from the floor

Before commencing procedure, introduce a pillow to protect the person's head from the hoist. Follow the procedures 1-2 as above.

Before undertaking stage 3, when using a mobile hoist, introduce the hoist by coming in from behind the head of the person where possible. Important - ensure brakes are on. Procedure 3 can now begin.

**To remove the sling, carry out the fitting procedure in reverse.
Important - always peel the sling away from the person to avoid skin friction and movement.**

Prism Hammock Sling - Material Variants / Uses

Standard Slings

Poly - The sling is made from a close knit polyester material and has nylon ripstop material to the outer side of the leg sections to aid the fitting of the sling. This is a general purpose sling and can be used for sitting and recumbent transfers between chairs, beds, baths, changing beds, W.C's etc.

Custom Built

Slipfit - The sling is made from a nylon ripstop (parachute style) material and has been designed to be easily fitted and removed. It is especially suited to situations where the client has limited ability to be moved to aid fitting the sling i.e. tight fitting chairs. This is a general purpose sling as the Poly but **cannot** be used in a bath as it will retain water around the client when being hoisted.

Spacer - The sling is made from knitted polyester "spacer" fabric that is manufactured with an 'air gap' internally. The material stretches slightly in all directions and moulds perfectly to the shape of the person being hoisted. This is a general purpose sling that be used as the Poly but has the added benefit of being suitable to be left in position if required (subject to risk assessment and monitoring of pressure risks).

Mesh - The sling is made from a polyester net material and has a non-absorbent foam in the leg sections. The inner surface of the leg sections has a close knit polyester covering for patient comfort. This sling is designed specifically for bathing and can be left in place whilst the person is being bathed / showered.

CE Marking - Declaration of Conformity

Registered Manufacturer No: CA003907

Prism UK Medical hereby declares that the Prism sling range meets the provisions of the European Union Council Directive 2007 / 47 / EEC concerning medical devices and are classified as Class 1 devices in accordance with that Directive.

Hammock Sling

Serial Number.....

We hereby certify that the products listed above are manufactured, labelled and packaged in accordance with the provisions of Directive 2007 / 47 / EEC and conform in all aspects to the specification(s) relevant to them at the time of despatch and are CE marked in accordance with Directive 2007 / 47 / EEC .



Signed : For Prism Medical UK

SLING CARE: LAUNDERING

GENERAL.

Sensitive laundering of lifting slings will extend the service life of the product. Conversely, improper laundering will reduce the service life and possibly affect the strength and integrity of the sling.

To be avoided:

- use of bleach at washing, particularly hypochlorite [may cause structural damage]
- overheating the sling [temperatures in excess of 100° may cause extreme shrinkage and structural damage]

Use of biological detergents is not recommended: though this type of detergent has no detrimental effect on the sling, the user's skin may be sensitive to any residue.

The use of fabric conditioners is not recommended.

Any hook & loop fastenings should be closed prior to laundering to avoid the 'hook' side damaging the sling or other items laundered with it.

GENERAL SOILING.

The sling should be laundered with soap or detergent at 50°C on a standard domestic cycle.

DISINFECTION.

The sling should be laundered to ensure a 'period at temperature' of:

- 10 minutes at 65°C
- 3 minutes at 72°C

Following NHS guidelines, extra time should be added to ensure thorough mixing of hot water through the wash load.

It may be that available equipment does not allow the required level of temperature control. The sling may be washed on a 90°C wash cycle, though it is likely that some moderate shrinkage and/or colour bleed will occur.

DRYING.

Slings will air dry quite quickly in a warm room.

Slings should not be placed over radiators or convectors.

Slings may be **CAREFULLY** tumble dried on a cool [60°C] setting and short cycle.

The sling should be checked at intervals not exceeding 5 minutes.

Do not allow the sling to overheat.

The sling should last in excess of 5 years if the above guidelines are adhered to. All Prism slings have a manufacturers guarantee of 12 months.

Care Plan Information

All slings in the Prism range have space on the labels to record the clients name. Each sling has unique individually colour coded lifting loops which helps maintain good lifting practice. The loop options facilitate variable positions from sitting to recline.

Date of Assessment :		FIM Score:	
Clients Name:			
Approx Weight:		D.O.B:	
Type of Sling:		Size:	
	Sling Straps Coloured loop Options		
Head Support (Option):	Adjust slider to suit the persons head position		
Shoulder Loop:	White <input type="checkbox"/> Orange <input type="checkbox"/> Lilac <input type="checkbox"/> Blue/black <input type="checkbox"/> Other <input type="checkbox"/>		
Leg loop:	Yellow <input type="checkbox"/> Green <input type="checkbox"/> Brown <input type="checkbox"/> Other <input type="checkbox"/>		
Navy inner leg strap does <u>NOT</u> attach direct to the carry bar at any time			
Comments			

If you have any questions regarding the use
of this equipment, please contact
Prism Medical UK, or your local authorised dealer.



Unit 1, Tir Llwyd Industrial Estate, St Asaph Avenue,
Kinmel Bay, Rhyl, Denbighshire, LL18 5JA

Tel +44 (0) 845 409 2000

www.prismmedical.co.uk

info@prismmedical.co.uk

This document conforms to BS EN ISO 10535:2006 Requirements